

ENDOMAX™

CORDLESS ENDODONTIC HANDPIECE



Order No.178

ISO 9001
ISO 13485



CE 0434

- Cordless convenience offers complete portability.
- Auto-stop and auto-reverse mechanisms for patient safety.
- Heavy-duty high torque motor.
- Reciprocating setting available without use of additional attachment.
- Ergonomic design provides comfort and easy operation.
- Features five torque settings.
- Automatic forward/reverse functions.
- Nine rotation speeds available for different applications.
- High capacity lithium-ion battery provides stable and prolonged usage.
- Offers quick and easy coupling.
- Low battery warning light.

DENTAMERICA®

ENDOMAX™

ENDOMAX™

CORDLESS ENDODONTIC HANDPIECE

Auto-Reverse Forward:



1. Clockwise rotation.
2. Reverses when torque limit is reached.
3. When the load is removed, rotation continues forward.

Reciprocating Mode:




The file reciprocates in a back and forth motion at 45°.

Auto-Reverse Stop:



1. Clockwise rotation.
2. Reverses when torque limit is reached.
3. When the load is removed, rotation stop.

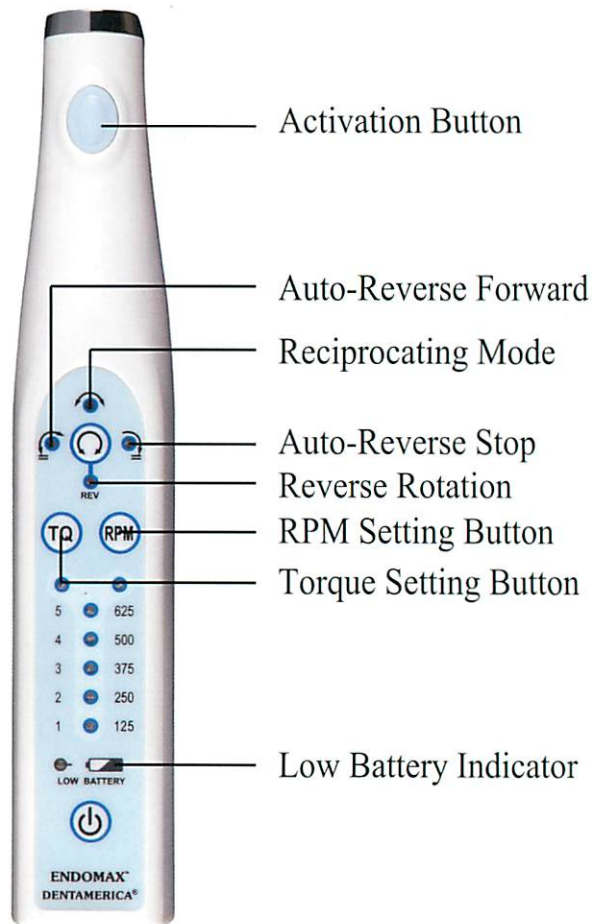

SPECIFICATIONS:

INPUT	AC 100V~240V 50/60Hz
FUSE	500mA 250V
BATTERY TYPE	3.7V 1620mAh
DIMENSIONS (HANDPIECE)	L192 xW26 xH29 mm
DIMENSIONS (CHARGING STATION)	L105xW85xH59mm
WEIGHT (HANDPIECE)	150g
WEIGHT (CHARGING STATION)	125g



Prestige Dental Products, Inc

3031 E. Coronado St., #D Anaheim, CA 92806

Phone # (877) 772-3888, Fax # (310) 388-8455

DENTAMERICA®