


 **apixia**[®]
DIGITAL IMAGING



apixia EXL[®] PSP SCANNER DIGITAL DENTAL RADIOGRAPHY IMAGING SYSTEM
minimum size. **maximum** performance.

FRESH IDEAS IN DIGITAL IMAGING

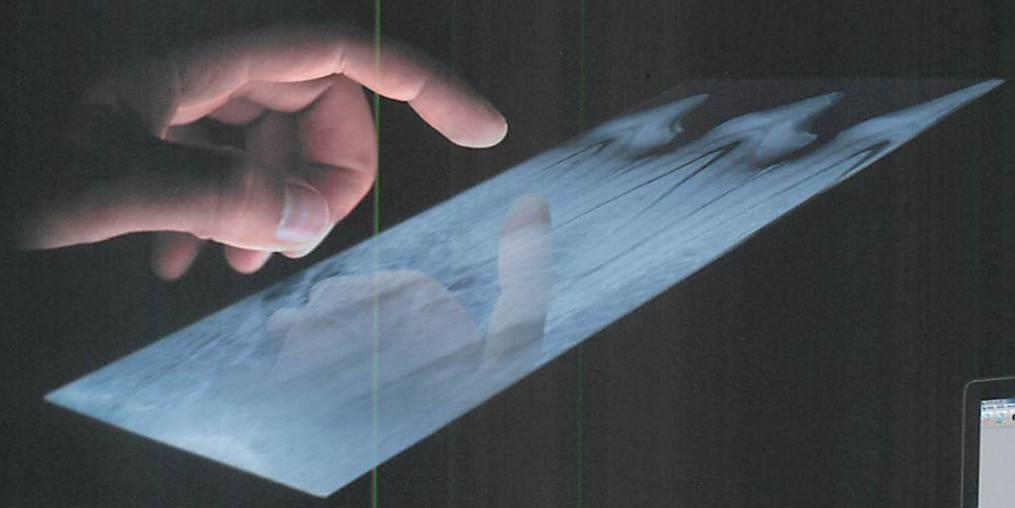


LASER FAST IMAGE PROCESSING

Imaging is fast and simple. The smart top-loading design harnesses the power of gravity for simple loading and fuss-free use. In just a short time, the crystal clear radiograph appears on the computer screen, ready for viewing. Simply insert the next exposed plate when ready.



INTELLIGENT

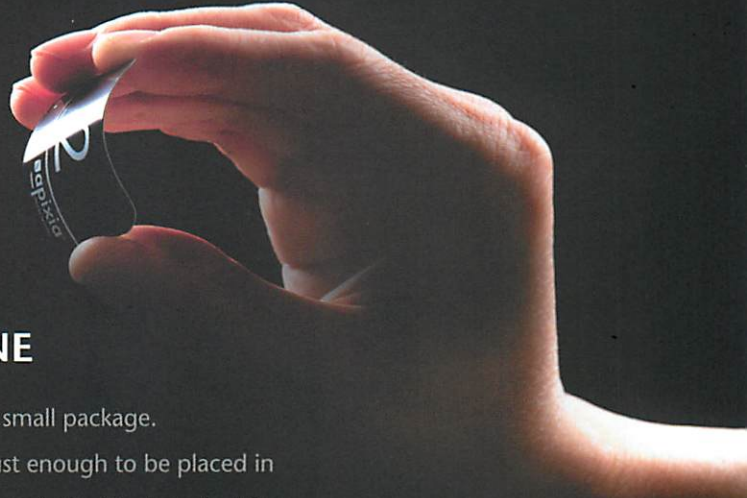

FLEXIBLE IN MORE WAYS THAN ONE

The apixia® PSP Scanner packs tremendous imaging power into a small package. Tiny enough to fit into the most crowded vet operatories, yet robust enough to be placed in a central location for use by the entire clinic, the apixia® PSP Scanner is the ideal product for any size vet practice. Imaging with PSP plates offers all the flexibility and durability of film, combined with the immediacy and environmental benefits of digital radiography.

REUSABLE

PSP plates are reusable a countless number of times to maximize your investment.

Every Size for Every Animal

apixia® PSP Scanner has been designed to take advantage of large plate sizes to accommodate larger size X-rays.




35x22mm



40x24mm



41x31mm

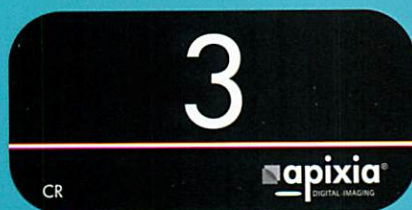

54x27mm



76x57mm

TECHNICAL SPECIFICATIONS

COMPUTER REQUIREMENTS*

Windows 10/Windows 8.1/Windows 8/Windows 7/
VISTA/XP (32, 64-bit)

USB 2.0 Port

CD/DVD Drive

*Computer not included

EXL® PSP SCANNER SPECIFICATIONS

Dimensions 214 mm x 170 mm x 266 mm

Weight 5.1 kg(11.23 lb)

Voltage AC 100-240V, 50-60Hz

Resolution 17 lp/mm

Scan Time 5-8 seconds

Connection Type USB 2.0 /TWAIN Compatible

Plate Sizes Size #0, #1, #2, #3, #4

Plate Erasing Automatic



Apixia Corporation Company Address: 7F-1, 128, Sec.2, Dunhua S. Rd., Da-An Dist., Taipei City Taiwan R.O.C.
Office: +886-2-27022608 Fax: +886-2-27554919 Email: sales@apixia.asia www.apixia.com

Apixia Inc. Company Address: 18688 San Jose Ave. Industry, CA 91748 USA
Office: 1-626-810-9201 Fax: 1-626-810-9207 Email: info@apixia.com www.apixia.net

Apixia Imaging Software offers optional integration to most leading brands of Practice Management Software with its patented DataGrabber®.

Contents Copyright 2015 Apixia Imaging, Inc. All Rights Reserved. Products are protected by numerous patents and pending patents worldwide. Apixia is independent (and thus not affiliated with or a licensee of) the Practice Management Software companies shown above. DataGrabber is a registered trademark of Apteryx, Inc. Product information subject to change without notice.